

**Bovine peroxisome proliferator activated receptor gamma (PPARG) ELISA Kit****Catalog No: EK760718**

<b>Reactivity</b>	:	Bovine
<b>Method Type</b>	:	Sandwich ELISA Detection
<b>Quantity</b>	:	96 tests
<b>Sample Type</b>	:	serum, plasma, Urine,tissue homogenates, cell culture supernates
<b>Detection Range</b>	:	0.31 ng/mL - 20 ng/mL
<b>Sensitivity</b>	:	0.31 ng/mL - 20 ng/mL

**Components:**

Assay plate (12 × 8 coated Microwells)	1
Standard: 22.5 ng/mL	1×0.5ml
Standard Diluent	1×1.5ml
HRP-Conjugate Reagent	1×6ml
Sample Diluent	1×6ml
Chromogen Solution A	1×6ml
Chromogen Solution B	1×6ml
Stop Solution	1×6ml
Wash Solution	1×20ml×30 fold
User manual	1
Adhesive Strip	1

**Product Principle:**

The kit is for the quantitative level of Bovine PPARG in the sample, adopt purified PPARG antibody to coat microtiter plate, make solid-phase antibody, then add PPARG to wells, Combine PPARG antibody with labeled HRP to form antibody-antigen - enzyme-antibody complex, after washing completely, add TMB substrate solution, TMB substrate becomes blue color at HRP enzyme-catalyzed, reaction is terminated by the addition of a stop solution and the color change is measured at a wavelength of 450 nm. The concentration of PPARG in the samples is then determined by comparing the O.D. of the samples to the standard curve.

**For research use only. Not for use in diagnostic procedures.**

### Specimen Requirements:

1. Serum-coagulation at room temperature for 10-20 min centrifuge at the speed of 2000-3000 rpm for 20-min. Remove supernatant, if precipitation appeared, Centrifuge again. Assay immediately or aliquot and store samples at -20°C or -80°C. Avoid repeated freeze-thaw cycles.
2. Plasma-use suited EDTA or citrate plasma as an anticoagulant, centrifuge at the speed of 2000-3000 rpm for 20-min. Remove supernatant, if precipitation appeared, centrifuge again. Assay immediately or aliquot and store samples at -20°C or -80°C. Avoid repeated freeze-thaw cycles. Centrifuge the sample again after thawing before the assay .
3. Urine-collect sue a sterile container, centrifuge at the speed of 2000-3000 rpm for 20-min. Remove supernatant, if precipitation appeared, Centrifuge again. The Operation of Hydrothorax and cerebrospinal fluid reference to it. Assay immediately or aliquot and store samples at -20°C or -80°C.
4. Cell culture supernatant-detect secretory components, Remove particulates by centrifugation for 20-min at the speed of 2000-3000 rpm. Remove supernatant detect the composition of cells, dilute cell suspension with PBS (PH7.2-7.4), Cell concentration reached 1 million / ml, repeated freeze-thaw cycles, damage cells and release of intracellular components, centrifugation 20-min at the speed of 2000-3000 rpm. remove supernatant, If precipitation appeared, Centrifugal again. Assay immediately or aliquot and store samples at -20°C or -80°C. Avoid repeated freeze-thaw cycles.
5. Tissue samples- After cutting samples, check the weight, Pipette PBS(PH7.2-7.4), Rapidly frozen with liquid nitrogen, maintainsamples at 2-8 after melting, Pipette PBS(PH7.4), homogenized by hand or Grinders, centrifugation 20-min at the speed of 2000-3000 rpm. Remove supernatant. Assay immediately or aliquot and store samples at -20°C or -80°C. Avoid repeated freeze-thaw cycles. Centrifuge the sample again after thawing before the assay.


Note:

1. Extract as soon as possible after Samples collection, and should be tested as soon as possible after the extraction. If not, samples must be stored at -20°C ( 1 month) or -80°C ( 2 months) to avoid loss of bioactivity and contamination.
2. Can' t detect the sample which contain NaN<sub>3</sub>, because NaN<sub>3</sub> inhibits HRP active.
2. Can' t detect the sample which contain NaN<sub>3</sub>, because NaN<sub>3</sub> inhibits HRP active.

### Reagent Preparation:

1: Wash Buffer (1x) - If crystals have formed in the concentrate, warm up to room temperature and mix gently until the crystals have completely dissolved. Dilute 20 ml of Wash Buffer Concentrate (30 x) into deionized or distilled water to prepare 600 ml of Wash Buffer (1 x).

2: Standard Dilute the standard Pipette 50l standard diluent in each tube. Pipette 100l standard (22.5 ng/mL) in the fifth tube. And take out 100l from the fifth five tube into the fourth. Pipette 50l from the fourth tube to the third tube and produce dilution series as below. The undiluted Standard serves as the high standard (22.5 ng/mL ). Sample Diluent serves as the zero standard blank well(0 ng/mL).

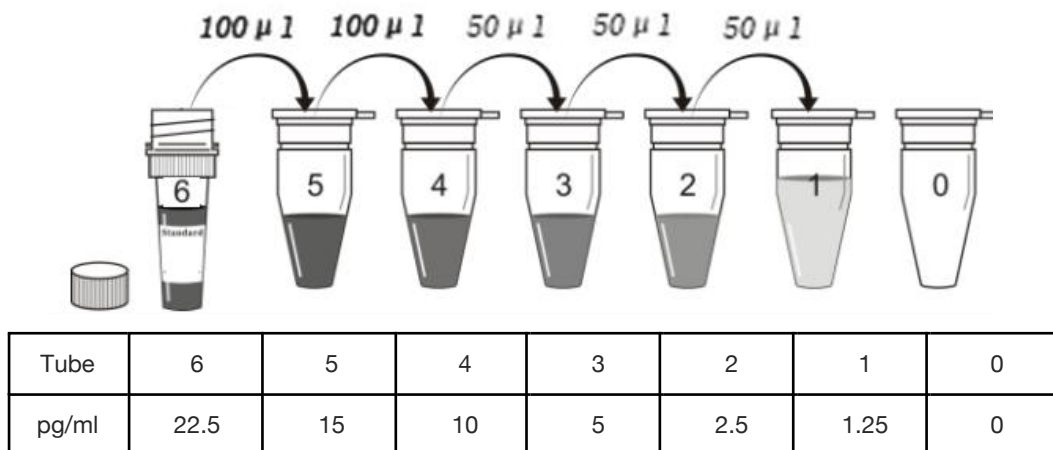

**For research use only. Not for use in diagnostic procedures.**

## Assay Procedure:

Step 1: Prepare all reagents, working standards, Blank and samples as directed in the previous sections.

Step 2: Refer to the Assay Layout Sheet to determine the number of wells to be used and put any remaining wells and the desiccant back into the pouch and seal the ziploc, store unused wells at 4°C

Step 3: Pipette standard 50l to testing standard well Pipette Sample diluent 40l to testing sample well, then add testing sample 10l (sample final dilution is 5-fold), Pipette sample to wells, don't touch the well wall as far as possible, and mix gently.

Step 4: Incubate: Cover with the adhesive strip provided, incubate for 30 min at 37.

Step 5: Configurate liquid: Dilute wash solution 30-fold with distilled water.

Step 6: Washing: Uncover the adhesive strip, discard liquid, pipette washing buffer to every well, still for 30s then drain, repeat 5 times.

Step 7: Add enzyme: Pipette HRP-Conjugate reagent 50l to each well, except blank well.

Step 8: Incubate: Operation with 4

Step 9: Washing: Operation with 6.

Step 10: Color: Pipette Chromogen Solution A 50ul and Chromogen Solution B 50ul to each well, avoid the light preservation for 15 min at 37.


Step 11: Stop the reaction: Pipette Stop Solution 50l to each well, stop the reaction (the blue change to yellow).

Step 12: Calculate: take blank well as zero. Read absorbance at 450nm after pipette Stop Solution within 15min.

## Calculation of Results:

Take the standard concentration as the horizontal, the OD value for the vertical, draw the standard curve on graph paper, Find out the corresponding concentration according to the sample OD value by the Sample curve, multiplied by the dilution multiple, or calculate the straight line regression equation of the standard curve with the standard concentration and the OD value, with the sample OD value in the equation, calculate the sample concentration, multiplied by the dilution factor, the result is the sample actual concentration.

Graphical Representation as following

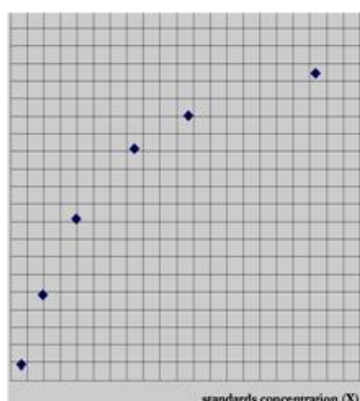

## Storage Conditions:

The unopened kit shall be stored at [2-8 ]

For opened kit can be stored at [2-8 ] for up to 1 month. If not be used recently, the standard should be kept in -20

**For research use only. Not for use in diagnostic procedures.**

**Precision:**

Intra-assay Precision (Precision within an assay): 3 samples with low, middle and high level Bovine PPARG were tested 20 times

on one plate, respectively.

Inter-assay Precision (Precision between assays): 3 samples with low, middle and high level Bovine PPARG were tested on 3 different plates, 8 replicates in each plate.

$CV(\%) = SD/mean \times 100$

Intra-Assay:  $CV < 8\%$

Inter-Assay:  $CV < 10\%$

**Attention:**

- 1: The kit takes out from the refrigeration should be balanced 15-30 minutes in the room temperature, if the coated ELISA plates have not been used up after opening, the plate should be stored in sealed bag.
- 2: Washing buffer will Crystallization separation, it can be heated in water to dissolve.
- 3: Pipette sample with pipettors each step, and proofread its accuracy frequently to avoid the experimental error. Pipette sample within 5 min, if the number of sample is big, recommend using multichannel pipettor.
- 4: If the testing material concentration is excessively high (The sample OD is higher than the first standard well) please dilute the sample (n-fold).
- 5: Adhesive Strip only limits the disposable use to avoid cross-contamination.
- 6: The substrate should evade the light to be preserved.
- 7: Please refer to the user instruction strictly, the test result determination must take the microtiter plate reader as a standard.
- 8: The preparation of samples and all the reagents should refer to infective material process.
- 9: Do not mix reagents with those from other lots.

**Washing Method:**

If you have technical questions about our products, please email us to [technical.support@afgsci.com](mailto:technical.support@afgsci.com) and have all your questions answered.